

Gran Torino


Biliary Endoprosthesis Kits



Hydrophilic polyurethane stent with 8 oval holes




Tapered pigtail tip with 4 oval holes

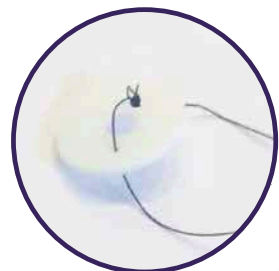

Button to be placed subcutaneously, against distal migration

GRAN TORINO Kit with hydrophilic polyurethane biliary endoprosthesis

12fr biliary stent, pigtail tip, for the long term internal drainage, with 13fr peel-away introducer system.

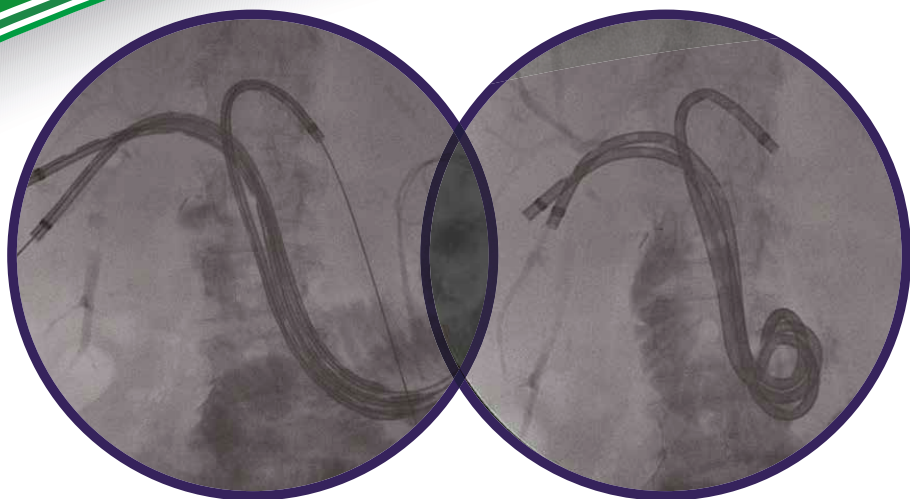

FEATURES

- Highly radiopaque, biocompatible, soft and encrustation resistant polyurethane
- Hydrophilic coating allows the easy introduction of the stent
- 8 large oval holes on the straight part and 4 on the pigtail, for the maximum drainage capacity
- Provided with an external catheter for coaxial use, to perform the initial drainage
- Proximal retention wire, connected to a special medical grade silicone subcutaneous button, to avoid the stent migration
- Kit contains a peel-away introducer and a pusher/straightener assembly

CE 1936

Cat. No. MB 856012 Gran Torino biliary endoprosthesis kit

DESCRIPTION	QUANT.
Biliary stent 12 Fr, 15 cm (straight portion) with button	1
Coaxial external drainage 6.6 Fr, 35cm, straight	1
Peel-away introducer system 13 Fr, 30/33 cm	1
Pusher/straightener	1

