



Hooker

Retrieving System

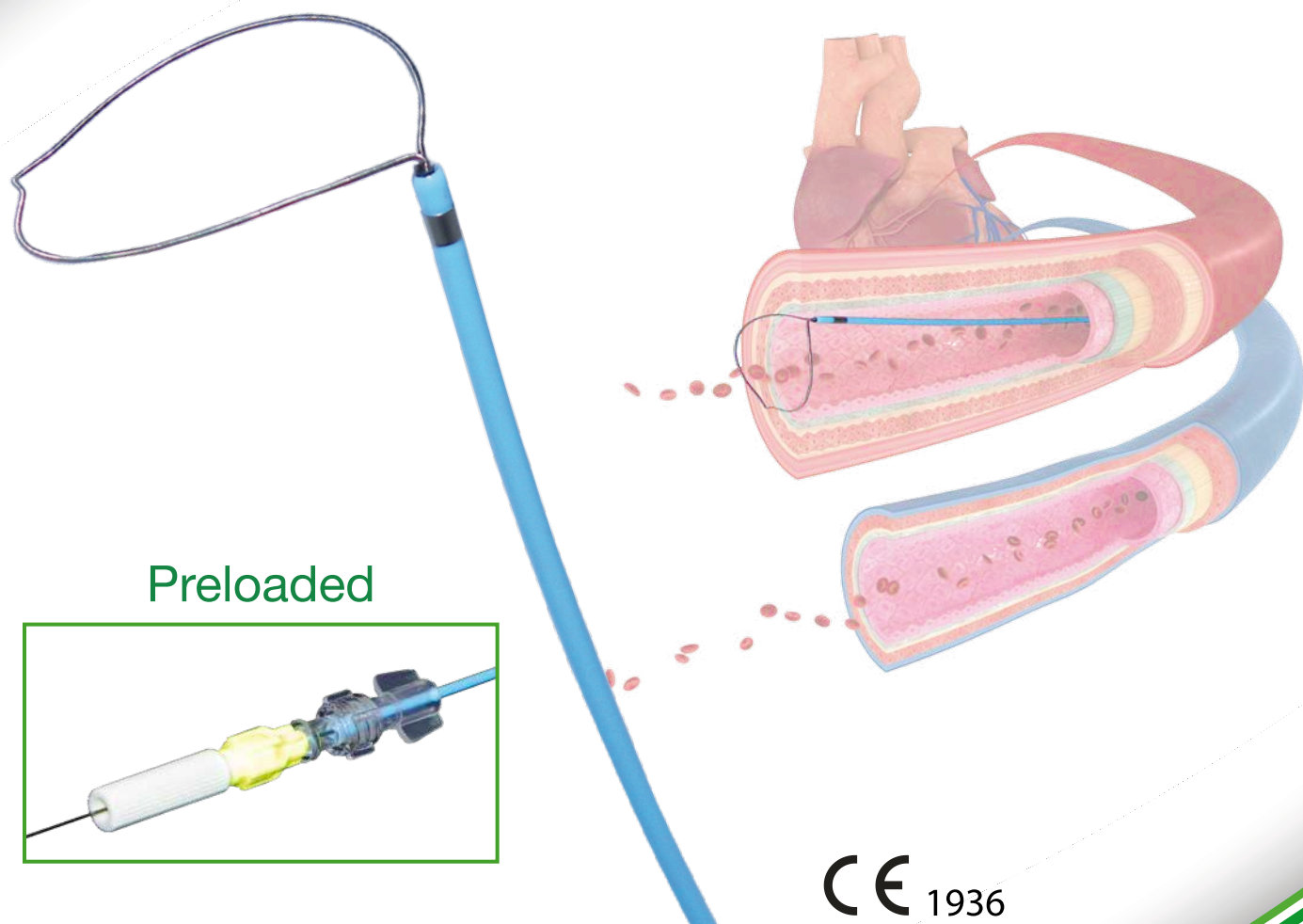

CE 1936

**Foreign body retrieval system
for vascular and non vascular
applications**

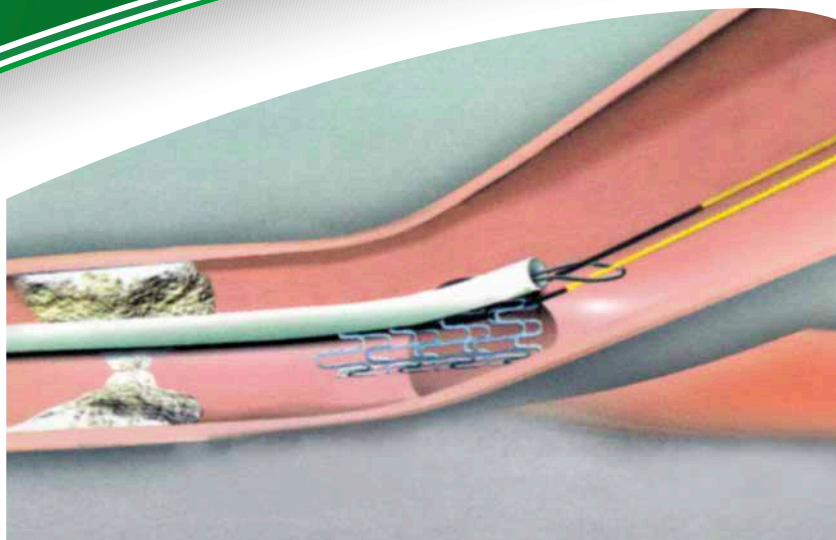

FEATURES

- The nitinol extruded loop provides superior resistance, flexibility, radiopacity and excellent shape memory, avoiding any deformation during the insertion procedure.
- The PTFE sheath allows the loop straightening and its passage through the catheter.
- The 90° loop to shaft enhances handling and provides an easy placement.
- Wide range of sizes for any requirements.
- Packaged with a guiding catheter to help the operator during the insertion step.
- Available as “microsnare” (loop diam. 3 and 6 mm), packaged with a guiding microcatheter and a “Y” connector.
- Sterile and disposable.

CODE NUMBER	CATHETER DIAMETER	LOOP DIAMETER	SNARE LENGTH
MB 807606	3.0 FR	3 mm	175 cm
MB 807608	3.0 FR	6 mm	175 cm
MB 807610	4.0 FR	6 mm	120 cm
MB 807612	5.0 FR	12 mm	120 cm
MB 807613	6.0 FR	18 mm	120 cm
MB 807614	6.0 FR	24 mm	120 cm
MB 807615	6.0 FR	30 mm	120 cm
MB 807616	6.0 FR	36 mm	120 cm
MB 807624	6.0 FR	24 mm	80 cm
MB 807625	6.0 FR	30 mm	220 cm

